

## BIOHIT COLONVIEW

quick test



### Product specifications

- Sample: stool sample
- Test principle: immunological detection of Hb and Hb/Hp complex
- Kit components: sampling tube, sampling paper, test cassette, instructions for use
- Storage: +2-30°C



In the USA and Japan for research use only.



### Ordering details:

Cat. No.	Item
602 250.02	30 tests/box

#### Biohit Oyj

Laippatie 1  
00880 Helsinki  
Finland  
Tel. +358 9 773 861  
Fax +358 9 773 86 200  
info@biohit.fi

#### Biohit HealthCare Srl

Via Boncompagni, 3  
20139 Milano  
Italy  
Tel +39 02 38238113  
Fax +39 02 38236521  
info.italy@biohit.fi

#### Biohit HealthCare Ltd.

Pioneer House, North Rd  
Ellesmere Port, CH65 1AD,  
United Kingdom  
Tel. +44 151 550 4 550  
Fax. +44 151 550 4 551  
info@biohithealthcare.co.uk

[www.biohithealthcare.com](http://www.biohithealthcare.com)

**BIOHIT HealthCare**  
Innovating for Health

**BIOHIT HealthCare**  
Innovating for Health



**BIOHIT COLONVIEW**  
quick test

VERSATILE SCREENING TOOL FOR COLORECTAL CANCER

490089EN 11/2014

RESULTS IN  
15 min

## BIOHIT COLONVIEW

quick test

### Detects bleeding derived from both lower and upper GI tract

Detection of occult blood in feces is essential in the early diagnosis of colorectal cancer. The Biohit ColonView quick test is an immunological lateral flow test which detects hemoglobin and hemoglobin-haptoglobin complex (degraded hemoglobin) from stool. With this combination of two biomarkers, the ColonView test can be used for detecting bleeding derived from the colon or upper parts of the gastrointestinal (GI) tract. This makes ColonView a superior tool to detect bleeding cancer tumors and their precursors.

- Detects bleeding from lower and upper GI tract
- Sensitivity 93-96 %, Specificity 80 %
- Hygienic sampling, transportation and analysis
- No dietary restrictions before sample collection
- Results in 15 min

### One sample for dual testing

Clever design of the ColonView sampling tube eliminates patient-errors and ensures that there will always be a correct sample amount. The sampling tube contains all the necessary reagents and no additional accessories are needed for the analysis. One sample is enough for the testing of both Hb and Hb/Hp complex, and the result is ready in 15 minutes.

### Testing made easy

The stool specimens are collected easily with the stool collection paper and sample applicator stick. To guarantee the best possible accuracy, it is recommended to take samples on 3 consecutive days.

1. Make bowel movement onto the collection paper


2. Unscrew the blue end and remove the sample applicator

3. Collect sample by inserting the stick into the stool specimen

4. Re-insert applicator stick

5. Unscrew the white cap and drop three drops to the test cassette

6. Read after 15 minutes

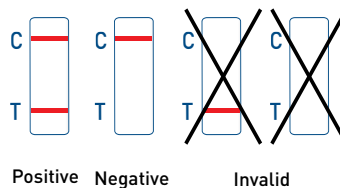

RESULTS IN  
15 min

### Improved accuracy to the screening of colorectal cancer

